

BAT STANDING-FRAME



BAT es un standing-frame que permite ir desde la posición horizontal hasta la posición vertical a través de mando eléctrico. Multiposición, permite assumir la posición prono o supino. Con sección de tronco y pélvis ajustable en altura a través de columna central, este standing-frame permite vários ajustes independientes al nivel de los soportes de rodillas y pies. Con función de abducción, desde 0° a 45° (totalizando 90°), ajustable de forma individual, con regulación de 5° em 5°. Soporte de cabeza 100% ajustable. Tablero ajustable en altura y inclinación. Cinturones de doble contención para mayor seguridad. Posibilidad de colocación de diversos accesorios.

BAT is a standing frame that allows to go from the horizontal to the vertical position through electric control. Multiposition, you can do prone or supine position. With trunk and pelvis section adjustable in height through central column, this standing frame also allows for several independent adjustments in knee and foot supports. With abduction function, from 0° to 45° (totalizing 90°) adjustable individually with regulation of 5° in 5°. 100% adjustable headrest. Board adjustable in height and inclination. Double containment belts for greater security. Several accessories available.

ORTHOSXXI, Unipessoal Lda.

Rua Santa Leocádia nº2735 // 4809-012 Sta. Leocádia de Briteiros // Guimarães // Portugal

Phone + 351 253 470 276 // Mobile +351 936 199 005 // Email margarida.silva@orthosxxi.com // Website www.orthosxxi.com

Características / Characteristics




Multiposición: prono y supino

Multiposition: prone and supine



Abducción hasta 45° (cada pierna)

Abduction until 45° (each leg)

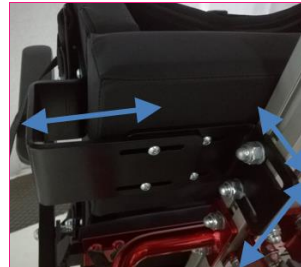

**Totalmente ajustable
Fully adjustable**



**Varios colores disponibles
Various colors available**

BAT standing-frame			
Marca/ Brand: Orthos XXI			
Longitud base Base length	1100mm	Soporte lateral tronco lateral trunk support	170x110mm
Ancho base Base width	600mm	Soporte lateral pélvis Pelvic lateral support	170x110mm
Soporte pélvico Pelvic support	380x180mm	Tablero (opcional) Board (optional)	650x550mm
Soporte de tronco Trunk support	380x180mm	Peso máximo usuario Maximum user weight	50Kg

Garantía/ Warranty
2 años contra defectos de fabricación 2 years against manufacturing defects
Tiempo estimado de entrega / Delivery time
Consultar/ Consult

ORTHOSXXI, Unipessoal Lda.

Rua Santa Leocádia nº2735 // 4809-012 Sta. Leocádia de Briteiros // Guimarães // Portugal

Phone + 351 253 470 276 // Mobile +351 936 199 005 // Email margarida.silva@orthosxxi.com // Website www.orthosxxi.com